


HEARTRAK SMART
INSTRUCTION SHEET

LIFE SUPPORT SYSTEMS
1-800-659-8151


GETTING STARTED

- * Placing the electrodes
The First electrode place on the upper right side of chest.
The Second electrode on the left lower side of chest.(See Diagram)
- * Connect the lead Wires.
(White snap to the upper right electrode, Black to the lower left)
- * Insert the lead-set cable firmly into the right side of the monitor.

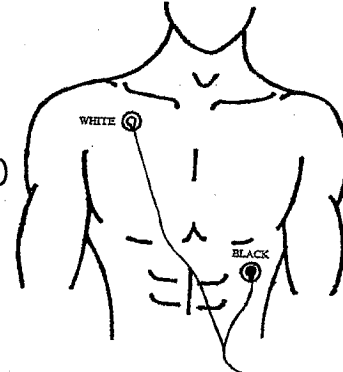


RECORDING AN EVENT/YOUR EKG

- * Press and Release the **RED RECORD** button found on the top of the monitor.
- * Remain quite and still. Breathe normally
- * Another event may be recorded if needed. Your Monitor can store _____ events.
(Please wait a full minute in between recordings)
- * If the monitor sends out a series of sequential beeps when trying to record an event
Please note the **Memory is full**.
- * You must call Life Support Systems first to clear out your recording.

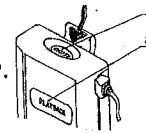


SENDING YOUR EKG OVER THE PHONE

- * Call Life Support Systems @ 1-800-659-8151
- * Slide the monitor out of the belt clip.(Pulls out of the top)
- * Upon direction from the staff:
- * Place the mouthpiece of your telephone over the small speaker hole found on the front of the monitor.
- * Monitor will produce a tone and is now playing back your stored recordings.
- * When the tone stops pick up the phone and await further instructions.

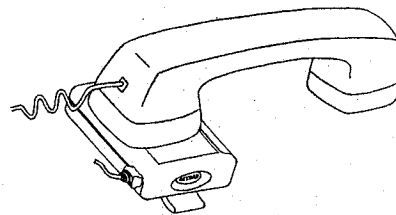

ERASING YOUR EVENTS

- * Press and hold the **PLAYBACK** button first while immediately pressing the **RECORD** button within two seconds a short beep is heard.
- * Your Monitor is now clear and ready for you to record more events.

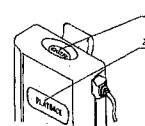

TO REPLACE BATTERIES

- * Remove the monitor out of the belt clip(Pulls out of the top).
- * Remove the sliding panel on the right side of the monitor.
- * Carefully remove the old battery. Replace with the new battery.
- * Battery is good for an average of 14 days.
- * **When the battery is low the monitor will produce an on-going series of sequential beeps.**

